




# AlloDerm™ is fully stocked.

We are ready to meet your needs no matter the size, big or small.

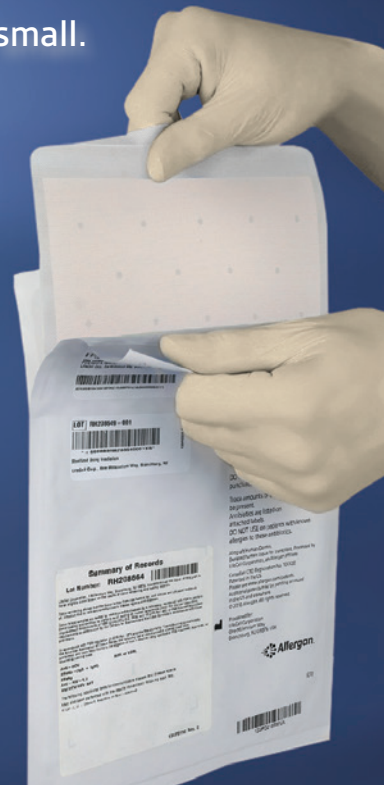

Order now by emailing [PRMOrder@allergan.com](mailto:PRMOrder@allergan.com) or by calling 1.800.367.5737.

To learn more, visit [AlloDerm.com](http://AlloDerm.com) and follow @AlloDermHCP

## INDICATIONS

ALLODERM SELECT™ Regenerative Tissue Matrix (ALLODERM SELECT™ RTM refers to both ALLODERM SELECT™ RTM and ALLODERM SELECT RESTORE™ RTM products) is intended to be used for repair or replacement of damaged or inadequate integumental tissue or for other homologous uses of human integument. This product is intended for single patient one-time use only. ALLODERM SELECT™ RTM is not indicated for use as a dural substitute or intended for use in veterinary applications.

## IMPORTANT SAFETY INFORMATION

### CONTRAINDICATIONS

ALLODERM SELECT™ RTM should not be used in patients with a known sensitivity to any of the antibiotics listed on the package and/or Polysorbate 20.

### WARNINGS

Processing of the tissue, laboratory testing, and careful donor screening minimize the risk of the donor tissue transmitting disease to the recipient patient. As with any processed donor tissue, ALLODERM SELECT™ RTM is not guaranteed to be free of all pathogens. No long-term studies have been conducted to evaluate the carcinogenic or mutagenic potential or reproductive impact of the clinical application of ALLODERM SELECT™ RTM.

**DO NOT** re-sterilize ALLODERM SELECT™ RTM. **DO NOT** reuse once the tissue graft has been removed from the packaging and/or is in contact with a patient. Discard all open and unused portions of the product in accordance with standard medical practice and institutional protocols for disposal of human tissue. Once a package or container seal has been compromised, the tissue shall be either transplanted, if appropriate, or otherwise discarded. **DO NOT** use if the foil pouch is opened or damaged. **DO NOT** use if the seal is broken or compromised. **DO NOT** use if the temperature monitoring device does not display "OK." **DO NOT** use after the expiration date noted on the label. Transfer ALLODERM SELECT™ RTM from the foil pouch aseptically. **DO NOT** place the foil pouch in the sterile field.

## PRECAUTIONS

Poor general medical condition or any pathology that would limit the blood supply and compromise healing should be considered when selecting patients for implanting ALLODERM SELECT™ RTM as such conditions may compromise successful clinical outcome. Whenever clinical circumstances require implantation in a site that is contaminated or infected, appropriate local and/or systemic anti-infective measures should be taken.

ALLODERM SELECT™ RTM has a distinct basement membrane (upper) and dermal surface (lower). When applied as an implant, it is recommended that the dermal side be placed against the most vascular tissue. Soak the tissue for a minimum of 2 minutes using a sterile basin and room temperature sterile saline or room temperature sterile lactated Ringer's solution to cover the tissue. If any hair is visible, remove using aseptic technique before implantation.

ALLODERM SELECT™ RTM should be hydrated and moist when the package is opened.

**DO NOT** use if this product is dry. Use of this product is limited to specific health professionals (e.g., physicians, dentists, and/or podiatrists). Certain considerations should be made to reduce the risk of adverse events when performing surgical procedures using a tissue graft. Please see the Instructions for Use (IFU) for more information on patient/product selection and surgical procedures involving tissue implantation before using ALLODERM SELECT™ RTM.

## ADVERSE EVENTS

The most commonly reported adverse events associated with the implant of a tissue graft include, but are not limited to the following: wound or systemic infection; seroma; dehiscence; hypersensitive, allergic or other immune response; and sloughing or failure of the graft.

**ALLODERM SELECT™ RTM is available by prescription only.**

For more information, please see the Instructions for Use (IFU) for ALLODERM SELECT™ RTM available at [www.allergan.com/AlloDermIFU](http://www.allergan.com/AlloDermIFU) or call 1.800.678.1605.

To report an adverse reaction, please call Allergan at 1.800.433.8871.

**Allergan**  
**Aesthetics**  
an AbbVie company

ALLODERM™, ALLODERM SELECT™, and their designs are trademarks of LifeCell Corporation, an AbbVie company.  
© 2021 AbbVie. All rights reserved. ALS144336 02/21



**Trust. Evidence. Experience. It's what AlloDerm™ RTM is made of.**

**No other ADM has been trusted with more than 2.5 million implantations.<sup>2</sup>**

**No other ADM has more publications, with hundreds of scientific\* and clinical articles.<sup>3</sup>**

**And no other ADM has the extensive experience of AlloDerm™ RTM, 25 years and counting.<sup>1</sup>**

\*Correlation of these results, based on animal studies, to results in humans has not been established.

#### INDICATIONS

ALLODERM SELECT™ Regenerative Tissue Matrix (ALLODERM SELECT™ RTM refers to both ALLODERM SELECT™ RTM and ALLODERM SELECT RESTORE™ RTM products) is intended to be used for repair or replacement of damaged or inadequate integumental tissue or for other homologous uses of human integument. This product is intended for single patient one-time use only. ALLODERM SELECT™ RTM is not indicated for use as a dural substitute or intended for use in veterinary applications.

#### IMPORTANT SAFETY INFORMATION

##### CONTRAINDICATIONS

ALLODERM SELECT™ RTM should not be used in patients with a known sensitivity to any of the antibiotics listed on the package and/or Polysorbate 20.

##### WARNINGS

**Processing of the tissue, laboratory testing, and careful donor screening minimize the risk of the donor tissue transmitting disease to the recipient patient. As with any processed donor tissue, ALLODERM SELECT™ RTM is not guaranteed to be free of all pathogens. No long-term studies have been conducted to evaluate the carcinogenic or mutagenic potential or reproductive impact of the clinical application of ALLODERM SELECT™ RTM.**

**DO NOT** re-sterilize ALLODERM SELECT™ RTM. **DO NOT** reuse once the tissue graft has been removed from the packaging and/or is in contact with a patient. Discard all open and unused portions of the product in accordance with standard medical practice and institutional protocols for disposal of human tissue. Once a package or container seal has been compromised, the tissue shall be either transplanted, if appropriate, or otherwise discarded. **DO NOT** use if the foil pouch is opened or damaged. **DO NOT** use if the seal is broken or compromised. **DO NOT** use if the temperature monitoring device does not display "OK." **DO NOT** use after the expiration date noted on the label. Transfer ALLODERM SELECT™ RTM from the foil pouch aseptically. **DO NOT** place the foil pouch in the sterile field.

##### PRECAUTIONS

Poor general medical condition or any pathology that would limit the blood supply and compromise healing should be considered when selecting patients for implanting

ALLODERM SELECT™ RTM as such conditions may compromise successful clinical outcome. Whenever clinical circumstances require implantation in a site that is contaminated or infected, appropriate local and/or systemic anti-infective measures should be taken.

ALLODERM SELECT™ RTM has a distinct basement membrane (upper) and dermal surface (lower). When applied as an implant, it is recommended that the dermal side be placed against the most vascular tissue. Soak the tissue for a minimum of 2 minutes using a sterile basin and room temperature sterile saline or room temperature sterile lactated Ringer's solution to cover the tissue. If any hair is visible, remove using aseptic technique before implantation.

ALLODERM SELECT™ RTM should be hydrated and moist when the package is opened. **DO NOT** use if this product is dry. Use of this product is limited to specific health professionals (e.g., physicians, dentists, and/or podiatrists). Certain considerations should be made to reduce the risk of adverse events when performing surgical procedures using a tissue graft. Please see the Instructions for Use (IFU) for more information on patient/product selection and surgical procedures involving tissue implantation before using ALLODERM SELECT™ RTM.

##### ADVERSE EVENTS

The most commonly reported adverse events associated with the implant of a tissue graft include, but are not limited to the following: wound or systemic infection; seroma; dehiscence; hypersensitive, allergic or other immune response; and sloughing or failure of the graft.

**ALLODERM SELECT™ RTM is available by prescription only.**

**For more information, please see the Instructions for Use (IFU) for ALLODERM SELECT™ RTM available at [www.allergan.com/AlloDermIFU](http://www.allergan.com/AlloDermIFU) or call 1.800.678.1605.**

**To report an adverse reaction, please call Allergan at 1.800.433.8871.**

**References:** 1. Wainwright DJ. Use of an acellular allograft dermal matrix (AlloDerm) in the management of full-thickness burns. *Burns*. 1995;21(4):243-248. 2. Data on file, Allergan. 2018. Sales Data. 3. Data on file, Allergan. PubMed search performed in June 2020.